

## Model E 28 | Drying and heating chambers with mechanical adjustment

With a BINDER E series heating oven, you obtain a reliable and powerful unit in a solid configuration at a very attractive price. This heating oven is preferred for tasks in basic research, as well as in human and veterinary medicine.

### BENEFITS

- Uniform drying conditions thanks to APT.line™ technology
- Identical test conditions throughout the chamber interior independent of sample size and quantity
- Outstanding thermal insulation saves operation costs



Model 28



Model 28

### IMPORTANT FEATURES

- Temperature range: +60 °C to +230 °C
- Adjustable exhaust air flap
- Hydro-mechanical thermostat
- Class 1 temperature limitation device
- Timer 0 – 120 min
- 2 chrome-plated racks

### ORDERING INFORMATION

Interior volume [L]	Power supply - unit fuse	Plug*	Version	Model version	Art.-No.
<b>Model E 28</b>					
28	230 V 1~ 50/60 Hz -	CEE 7/7	Standard	E028-230V	9010-0001

## TECHNICAL DATA


Data	
Article Number	9010-0001
Designation	E028-230V
Option model	Standard
Performance Data Temperature	
Temperature range	60...230 °C
Heating up time to 150°C	36 min
Recovery time after door was opened for 30 s at 150°C	19 min
Electrical data	
Rated Voltage	230 V
Power frequency	50/60 Hz
Nominal power	0,8 kW
Phase (Nominal voltage)	1~
Dimensions and weights	
Interior volume	28 L
Net weight of the unit (empty)	22 kg
Load per rack	10 kg
Permitted load	25 kg
Wall clearance back	160 mm
Wall clearance sidewise	100 mm
Housing dimensions not incl. fittings and connections	
Width net	580 mm
Height net	405 mm
Depth net	425 mm
Internal Dimensions	
Interior width	400 mm
Interior height	280 mm
Interior depth	250 mm
Unit doors	1
Fixtures	

Number of shelves (std./max.)

2/4

All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10%. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

## DIMENSIONS Incl. fittings and connections [mm]

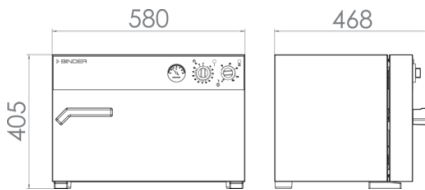

## OPTIONS

Designation	Description	E 28	*	Art.-No.
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	•	-	8012-1114
Calibration certificate, temperature	for temperature, measurement in center of chamber at specified temperature	•	-	8012-1133

## ACCESSORIES

Designation	Description	E 28	*	Art.-No.
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	-	8012-2250
Rack	chrome plated	•	-	8012-2032
Rubber pads	set anti-slip feet	•	-	8012-0702
Shelf, perforated	Stainless steel	•	-	8012-2160

## SERVICES

Designation	Description	*	Art.-No.
<b>Calibration services</b>			
Calibration certificate, temperature	Calibration of one (1) test temperature specified by the user in center of chamber, including certificate	05	DL30-0110
	Extension of calibration of one (1) additional test temperature specified by the user in the center of the usable space, including certificate	05	DL30-0102
<b>Installation services</b>			
Unit commissioning	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work)	05	DL10-0160
Unit instructions	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming	05	DL10-0530
<b>Maintenance contracts</b>			

Designation	Description	*	Art.-No.
BRONZE 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts	05	DL20-0720
GOLD 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, replacement of wear parts, calibration of one test temperature specified by the user in the center of the usable space, including certificate	05	DL20-0915
SILVER 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, calibration of one test temperature specified by the user in the center of the usable space, without certificate	05	DL20-0815
<b>Maintenance services</b>			
Maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Calibration of a test temperature specified by the user in center of usable space without certificate	05	DL20-0601
<b>Warranty service</b>			
1-year warranty extension	The warranty is extended by 1 year from the delivery date, wear parts are excluded	05	DL50-0010

**NOTES**

- 05 Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

**BINDER GmbH**

Tuttlingen, Germany  
TEL +49 7462 2005 0  
info@binder-world.com  
www.binder-world.com

**BINDER Inc.**

Bohemia, NY, USA  
TEL +1 631 224 4340  
usa@binder-world.com  
www.binder-world.us

**BINDER Environmental Testing  
Equipment (Shanghai) Co., Ltd.**

Shanghai, P.R. China  
TEL +86 21 685 808 25  
china@binder-world.com  
www.binder-world.com