


## Model FD-S 56 | Drying and heating chambers with forced convection

The new products in the Solid.Line impress with their tried-and-tested BINDER quality and reliability. They are suited to drying and heating. The attributes of the Solid.Line are geared towards the application in question and therefore used in research and quality assurance.

### BENEFITS

- Easy handling
- Accurate temperature control
- Simple and ergonomic door opening

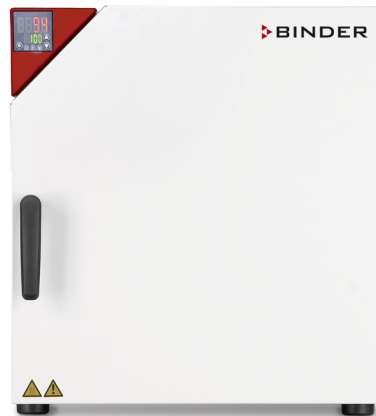

Model 56

### IMPORTANT FEATURES

- Temperature range: +10°C above ambient temperature to +250 °C
- APT.line™ preheating chamber technology
- Forced convection
- Adjustable exhaust air flap
- Controller with timer function
- 1 chrome-plated rack, incl. shelf supports
- Class 2 integrated independent adjustable temperature safety device (DIN 12880) with visual alarm

### ORDERING INFORMATION


Interior volume [L]	Power supply - unit fuse	Plug*	Version	Model version	Art.-No.
<b>Model FD-S 56</b>					
55	230 V 1~ 50/60 Hz -6,3 A	CEE 7/7	Standard	FDS056-230V	9090-0018
	120 V 1~ 60 Hz -6,3 A	NEMA 5-15	Standard	FDS056UL-120V	9090-0019

## TECHNICAL DATA

Data		
Article Number	9090-0018	9090-0019
Designation	FDS056-230V	FDS056UL-120V
Option model	Standard	Standard
Performance Data Temperature		
Temperature range	+10 °C above ambient temperature to 250 °C	+10 °C above ambient temperature to 250 °C
Temperature uniformity at 150°C	2.6 ± K	2.6 ± K
Temperature fluctuation at 150°C	0.4 ± K	0.4 ± K
Heating up time to 150°C	15 min	15 min
Recovery time after door was opened for 30 s at 150°C	5 min	5 min
Air change data		
Air change rate at 100°C	4 x/h	4 x/h
Electrical data		
Rated Voltage	230 V	120 V
Power frequency	50/60 Hz	60 Hz
Nominal power	1,1 kW	1,2 kW
Unit fuse	6,3 A	6,3 A
Phase (Nominal voltage)	1~	1~
Dimensions and weights		
Interior volume	55 L	55 L
Net weight of the unit (empty)	35 kg	35 kg
Load per rack	15 kg	15 kg
Permitted load	45 kg	45 kg
Wall clearance back	160 mm	160 mm
Wall clearance sidewise	100 mm	100 mm
Housing dimensions not incl. fittings and connections		
Width net	614 mm	614 mm
Height net	625 mm	625 mm
Depth net	577 mm	577 mm
Internal Dimensions		

Interior width	400 mm	400 mm
Interior height	440 mm	440 mm
Interior depth	310 mm	310 mm
Unit doors	1	1
<b>Environment-specific data</b>		
Sound-pressure level	43 dB(A)	43 dB(A)
Energy consumption at 150°C	290 Wh/h	290 Wh/h
<b>Fixtures</b>		
Number of shelves (std./max.)	1/3	1/3

### DIMENSIONS Incl. fittings and connections [mm]

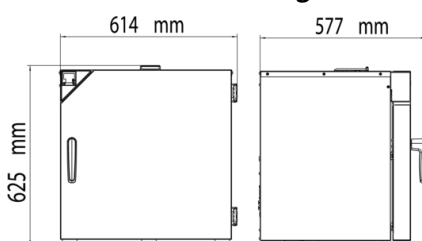

### ACCESSORIES

Designation	Description	FD-S 56	*	Art.-No.
Grid shelf	Incl. shelf supports			
	chrome plated	•	-	8012-1807
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	-	8012-2250
Rubber pads	set anti-slip feet	•	-	8012-2030

### SERVICES

Designation	Description	*	Art.-No.
<b>Calibration services</b>			
Calibration certificate, temperature	Calibration of one (1) test temperature specified by the user in center of chamber, including certificate	05	DL30-0120
	Extension of calibration of one (1) additional test temperature specified by the user in the center of the usable space, including certificate	05	DL30-0102
Measurement of air exchange rate	including certificate (in accordance with ASTM D5374)	05	DL33-0000
Temperature measurement, 9 measuring points	Temperature measurement with 9 measuring points with a set value specified by the user, including certificate	05	DL30-0109

Designation	Description	*	Art.-No.
Temperature measurement, 27 measuring points	Temperature measurement with 27 measuring points with a set value specified by the user, including certificate	05	DL30-0127
<b>Installation services</b>			
Unit commissioning	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work)	05	DL10-0160
Unit instructions	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming	05	DL10-0530
<b>Maintenance contracts</b>			
BRONZE 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts	05	DL20-0720
GOLD 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, replacement of wear parts, calibration of one test temperature specified by the user in the center of the usable space, including certificate	05	DL20-0910
SILVER 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, calibration of one test temperature specified by the user in the center of the usable space, without certificate	05	DL20-0810
<b>Maintenance services</b>			
Maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Calibration of a test temperature specified by the user in center of usable space without certificate	05	DL20-0602
<b>Warranty service</b>			
1-year warranty extension	The warranty is extended by 1 year from the delivery date, wear parts are excluded	05	DL50-0010

**NOTES**

- 05 Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

**BINDER GmbH**

Tuttlingen, Germany  
TEL +49 7462 2005 0  
info@binder-world.com  
www.binder-world.com

**BINDER Inc.**

Bohemia, NY, USA  
TEL +1 631 224 4340  
usa@binder-world.com  
www.binder-world.us

**BINDER Environmental Testing  
Equipment (Shanghai) Co., Ltd.**

Shanghai, P.R. China  
TEL +86 21 685 808 25  
china@binder-world.com  
www.binder-world.com